



Specifications

Power Supply	100~240VAC, 50/60Hz, 500VA
Flow Rate	0.1~45ml/s
Pressure Limit	1200PSI
Auto Filling Rate	8.0ml/s
Air Purge Rate	0.5~9.9ml/s
Piston Rod Speed	9.9ml/s
Rise Time	0.1-9.9sec
Injection Delay	0.0-99.9sec
X-Ray Delay	0.0-99.9sec
Phases	4
Injection Programs	500
Injection Records	500
Volume	1~150ml
Multiple Language User Interface	

High Quality Disposables



· 150ml Syringe for Contrast



ANT000150

Dual Display Operation System



DSA Injector

Leading Contrast Media Injectors Forward

antmed

Shenzhen Antmed Co., Ltd.

18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China
 Tel: +86 755 8606 0992 Fax: +86 755 8606 0991
 E-mail: info@antmed.com http://www.antmed.com



P/N: EN-DSA-201912

antmed
 Partner For Life

DSA Injector

Introducing an all-new ImaStar DSA Injector, offering an extended list of product features, reliable performance, and a safe patient experience.



Filling interface & parameter setting interface intelligent switching

Cantilever

- Large cantilever allows for extended reach to ensure smooth workflow
- 28" (72 cm) swing radius and 11" (28 cm) vertical height adjustment range

Intelligent Design

- Automatic head angle detection to prevent false start
- Touch-screen on the injector can toggle between filling interface and injecting interface automatically. When the injector head is upward, the filling interface is displayed when the injector head is downward.
- Air purge confirmation prompt is on the injector head touch-screen display

Heat Maintainer

- Ensures contrast maintains is at ideal temperature

Connect to Other Imaging Equipment

- Connect with other imaging equipment for enhanced imaging use contrast more efficiently and minimize radiation

Color Touch-Screen

- 12.1" and 12.1" color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence

Injection Record Storage


- Stores up to 500 injection history records



Bluetooth communication

Injection Ready Status

- Choose between single injection or multiple injection program

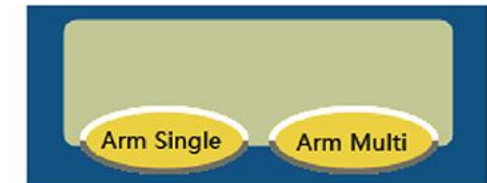

A Variety of Configuration Options

- Stores up to 500 customize injection programs
- Program storage capacity can be increased based on user needs



Multiple Operation Functions

- Offers 4 injection phases to meet demands of complex imaging needs

Safety

- To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- Automatically slows down the injection when approaching to pressure limit
- Injection stops when pressure limit is reached. Alerts with both light and sound
- Waterproof design inside the head to minimize damage from contrast media and saline leakage