

Specifications

Power Supply	100~240VAC, 50/60Hz, 200VA
Flow Rate	0.1~10ml/s
Pressure Limit	325PSI
Auto Filling Rate	8ml/s
Air Purge Rate	1.0~9.9ml/s
Piston Rod Speed	9.9ml/s
Phases	6
Injection Programs	2000
Injection Records	2000
Volume	A: 1~200ml
	B: 1~200ml
Multiple Language Us	er Interface

High Quality Disposables



- · 200ml Syringe for Contrast
- · 200ml Syringe for Saline

antmed"

Shenzhen Antmed Co., Ltd. 18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China

Tel: +86 755 8606 0992 Fax: +86 755 8606 0991 E-mail: info@antmed.com http://www.antmed.com






CT Dual Injector

Leading Contrast Media Injectors Forward



CT Dual Injector

Introducing an all-new ImaStar CT dual injector, offering an extended list of product features, reliable performance, and a safe patient experience.

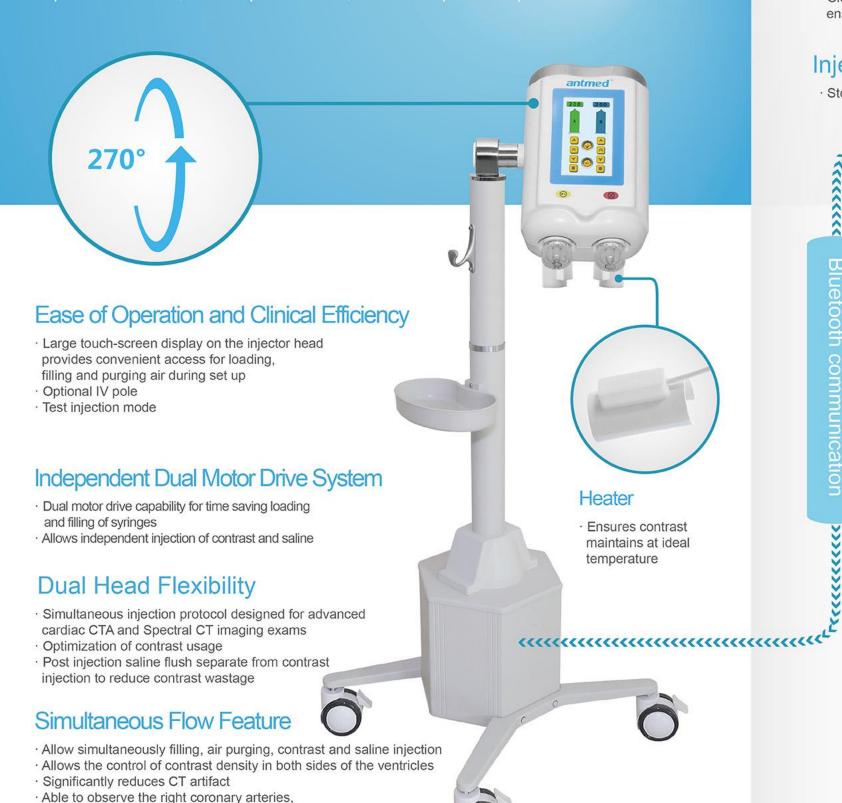

right atrium and right ventricle in the same image

Color Touch-Screen

- · 12.1 and 7" color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence

Injection Record Storage



Auto Detect

· Syringe detection feature automates air purging and piston rod operation

Connect to Other Imaging Equipment

- · Connect with other imaging equipment for enhanced imaging, use contrast more efficiently and minimize radiation
- Option

Bluetooth communication


Heater

Ensures contrast

maintains at ideal temperature

Real-time Pressure Curve

· Display provide real-time pressure monitoring curve

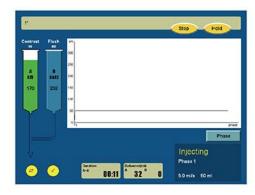

A Variety of **Configuration Options**

- · Stores up to 2000 customized injection
- Program storage capacity can be increased based on user needs



Safety

- · To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- · Automatically slows down the injection when approaching to pressure limit
- · Injection stops when pressure limit is reached. Alerts with both light and sound
- · Waterproof injector head designed to minimize damage from contrast and saline leakage
- · Ability to pause and resume while injection is underway